

SOMNOtouch™ PSG

OUR BIGGEST IDEAS. NOW IN OUR SMALLEST PSG DEVICE

Headbox with Continuous Impedance Check

Lower noise, higher signal quality and easy confirmation of proper electrode hook-up

AASM Compliant Recording (up to 25 Channels)

Headbox with EEG, EOG, EMG. Combi sensor for PLM, ECG & effort. 11 Internal channels inc: SpO₂, pulse rate & plethysmogram

Real-time Telemetric Data Transfer

On a PC, Smartphone or Tablet

55 mm

Cuffless Continuous Blood Pressure Measurement

A new paradigm in measuring blood pressure in a sleep lab - without disturbing sleep (optional upgrade)

84 mm

Intelligent Connect™

Automatic sensor recognition - making montages a thing of the past

Small Compact Robust Design


Credit card footprint means unparalleled comfort. Touch Screen display for PC independent set-up & signal quality check at bedside

ACTUAL
SIZE

The **SOMNOtouch™ PSG** combines the comfort, function and power of much larger sleep diagnostic devices in a small perfectly formed size. At only 64 g and 84 mm x 55 mm it is one of the lightest yet powerful touch screen ambulatory/stationary polysomnography devices on the market today.

Colour Touch Screen

- ◆ Start the data recording on the device
- ◆ Signal check on the screen
- ◆ Zoom to a selectable time base
- ◆ Easy to use interface
- ◆ Clear icons to easily identify attached sensors

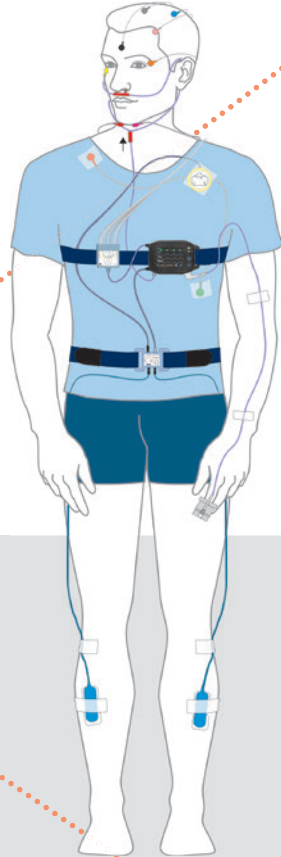


AUTOMATIC SOFTWARE & FIRMWARE UPDATES



Versatility – In a compact design

- ◆ Respiratory screening, plus recording of EEG, EMG & EOG signals from the headbox and of the ECG & PLM signals with a combi-sensor
- ◆ Minimal costs per measurement
- ◆ Configurable across reimbursement systems

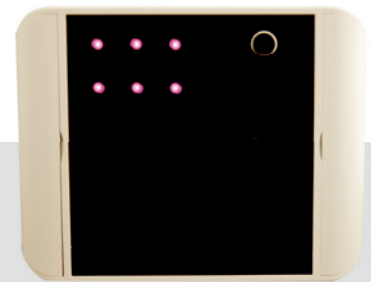

Compact Camera, Item no.: BTV101

Intelligent Connect™

- ◆ Automatic sensor recognition, no montages required
- ◆ Clear icons to easily identify functions/sensors
- ◆ AASM compliant headbox with continuous impedance recording

Video & Accessories (optional)

- ◆ Now with on-line recording capabilities
- ◆ Powered by our powerful DOMINO Software
- ◆ Complete range of stationary video solutions
- ◆ Comfortable Bluetooth receiver (BTR062)
- ◆ Include external devices via our ADC board

Never miss a new feature or update

Docking Station

- ◆ Fast recharging for PSG recording
- ◆ Fast data transfer (Approx. 90s for a PSG)
- ◆ Virtual dockingstation to display device status on PC
- ◆ Available in a 3 port design for bigger installations

